



Sample Name : M3194-Nocodazole-09
Acq. Operator : C22
Spec. Reported : MS Integration
Acq. Method : D:\METHODS\1-POS-3MIN.M
Analysis Method : D:\METHODS\1-POS-3MIN.M
Sample Info : Easy-Access Method: '1-POS-3MIN.M'
Method Info : Column:Xbridge C18(2) (4.6x 50 mm, 3.5um)
Mobile phase: H2O(10 mmol NH4HCO3) (A) / ACN(B)
Elution program: Gradient from 5 to 95% of B in1.5min at 1.8ml/min
Temperature: 50 °C
Detection: UV (214 , 4 nm) and MS (ESI, Pos mode ,103 to 800 amu)


*DAD1 A, Sig=214,4 Ref=off

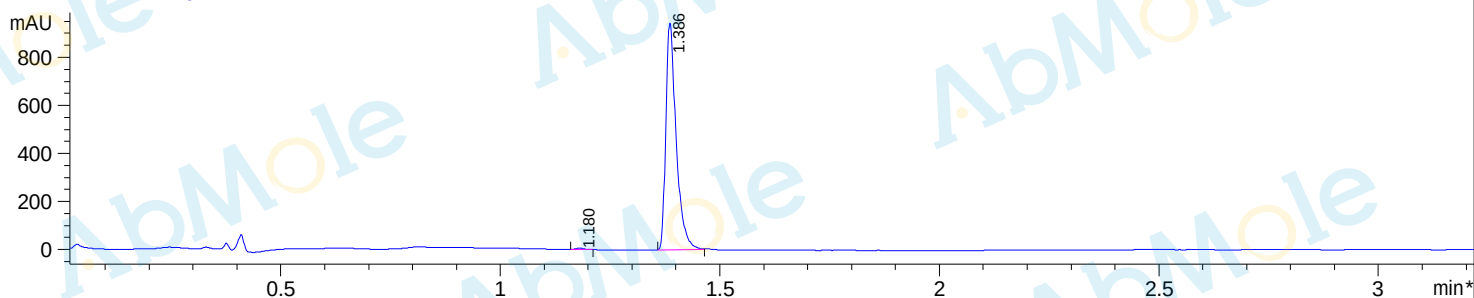

*DAD1 B, Sig=254,4 Ref=off

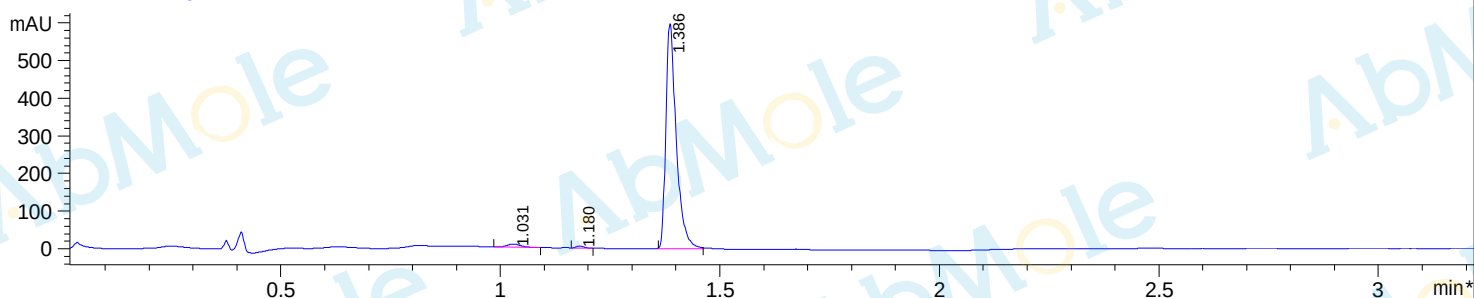

MSD1 TIC, MS File ES-API, Pos, Scan, Frag: 80

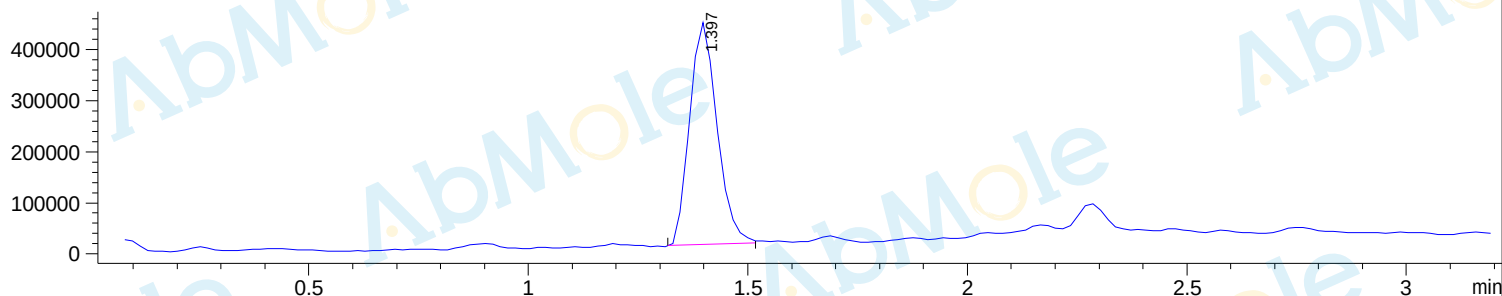

Integration Results for DAD1 A, Sig=214,4 Ref=off

RetTim	Width	Area	Height	Area%
1.18	0.02	7.97	6.39	0.51
1.39	0.02	1554.74	945.09	99.49

Integration Results for DAD1 B, Sig=254,4 Ref=off

RetTim	Width	Area	Height	Area%
1.03	0.03	17.56	8.21	1.74
1.18	0.02	6.73	5.47	0.67
1.39	0.02	985.35	598.48	97.59

Ret. Time: 1.40

<<<< POSITIVE SPECTRA >>>>

*ES-API Positive

