

LCMS REPORT

Sample Information


LCMS

Acquired by : System Administrator Date Acquired : 2021/5/15 14:13:25
Sample Name : Evans Blue-01 injection Volume : 3
Method File : B 5%-95%-P-new - 2.6min.lcm Vial : 9


Column : Waters SunFire C18 50*4.6mm 5um 2.000ml/min 2.6min
Pump A : 0.1% FA and 10% ACN in H2O
Pump B : 0.1% FA and 10% H2O in ACN

Chromatogram


mAU

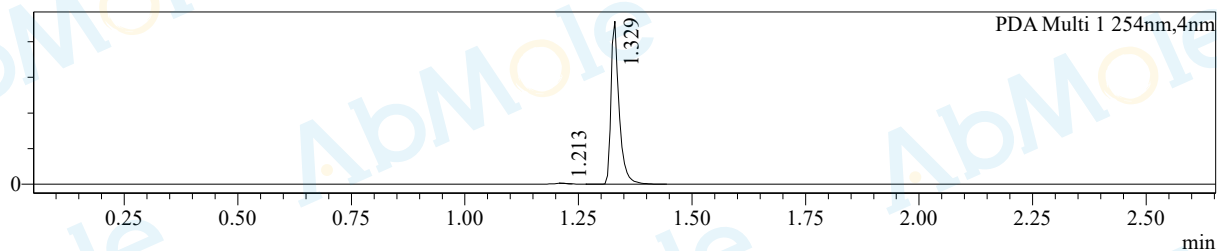

mAU

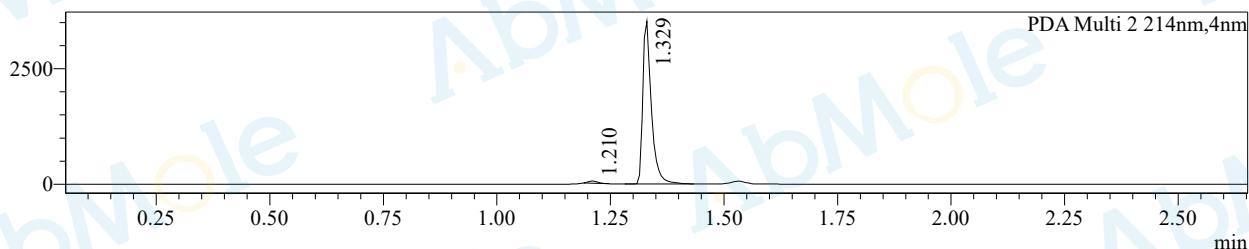

段号1

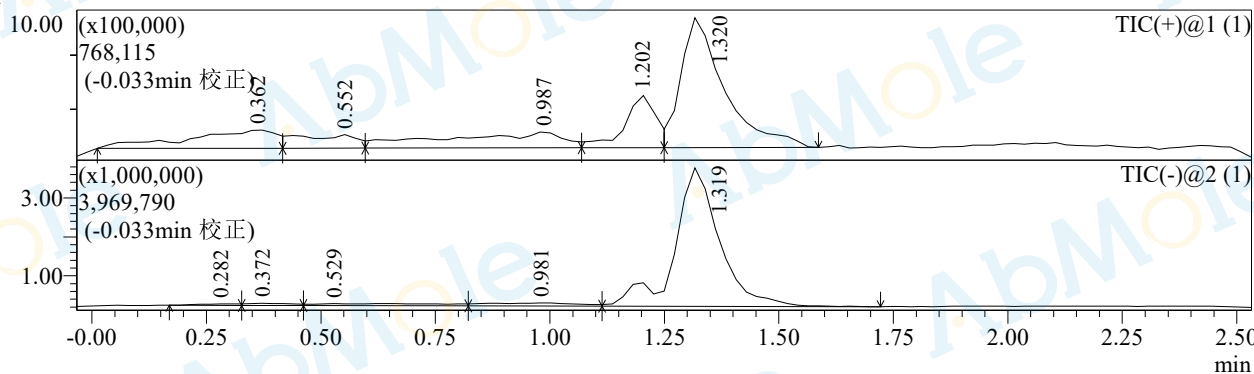

Integration Results for DAD

PDA Ch1 254nm


Peak#	Ret. Time	Height	Area	Area%
1	1.213	3703	2120	0.074
2	1.329	2290022	2857440	99.926
总计		2293725	2859559	100.000

PDA Ch2 214nm


Peak#	Ret. Time	Height	Area	Area%
1	1.210	44145	67887	1.430
2	1.329	3526216	4679099	98.570
总计		3570361	4746987	100.000

Mass Spectrum


Rt=0.287min 负

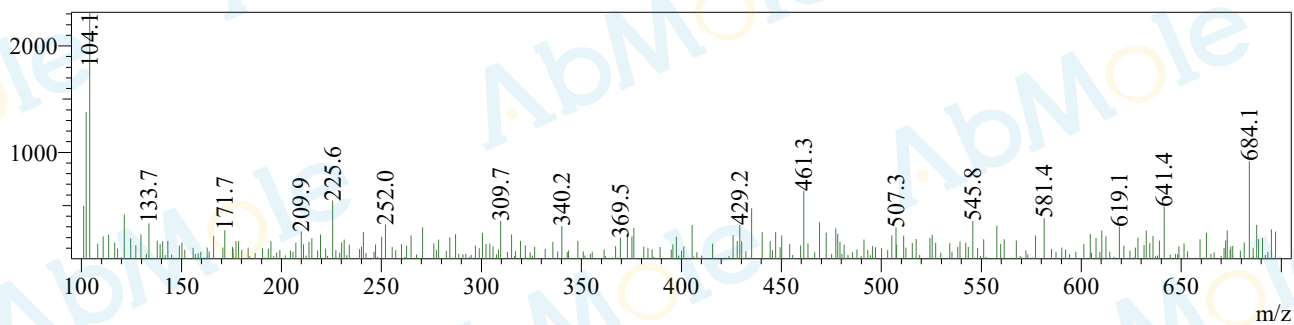


Rt=0.349min 正

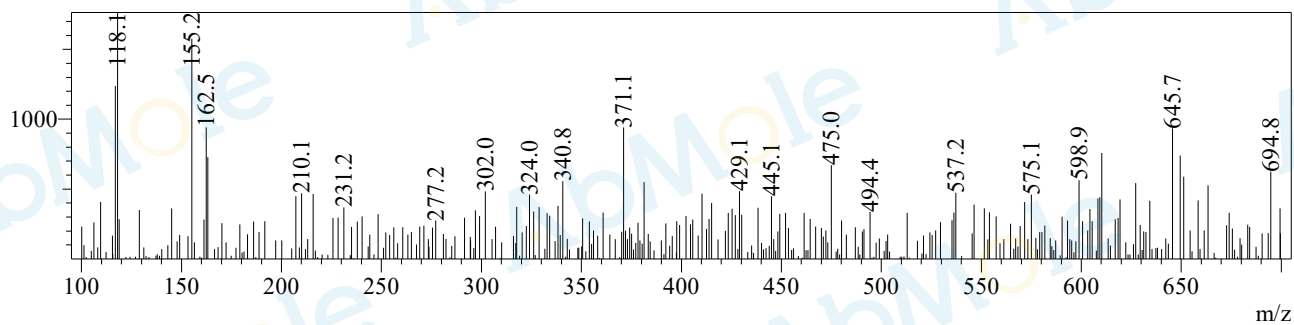


Rt=0.377min 负

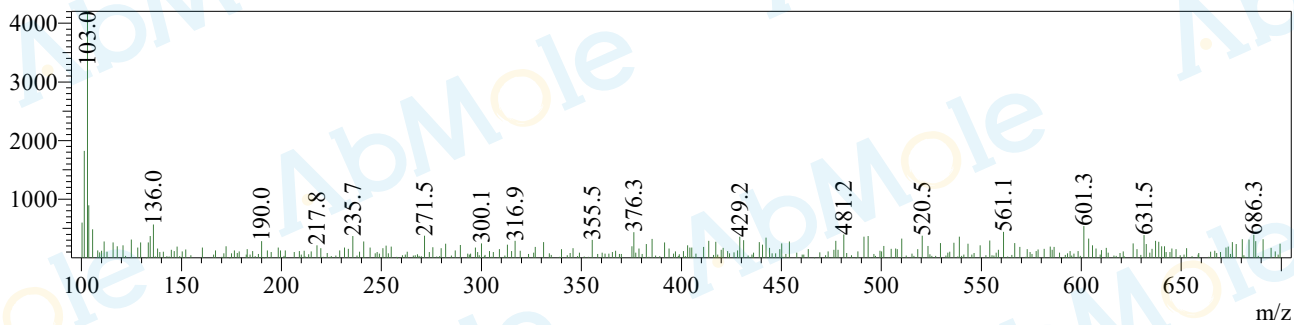


Rt=0.535min 负

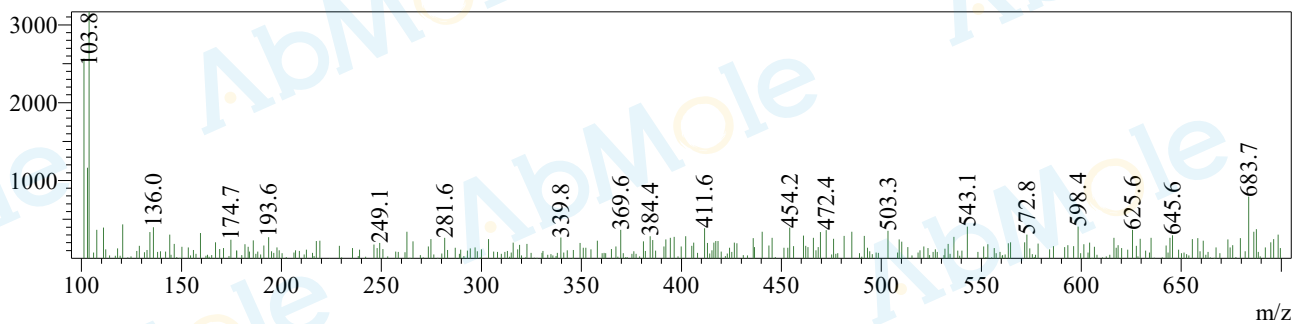


Rt=0.552min 正

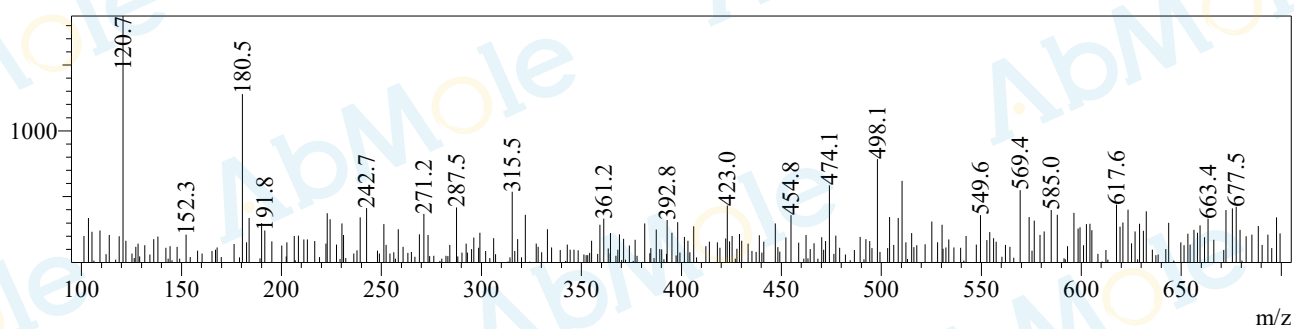


Rt=0.985min 负

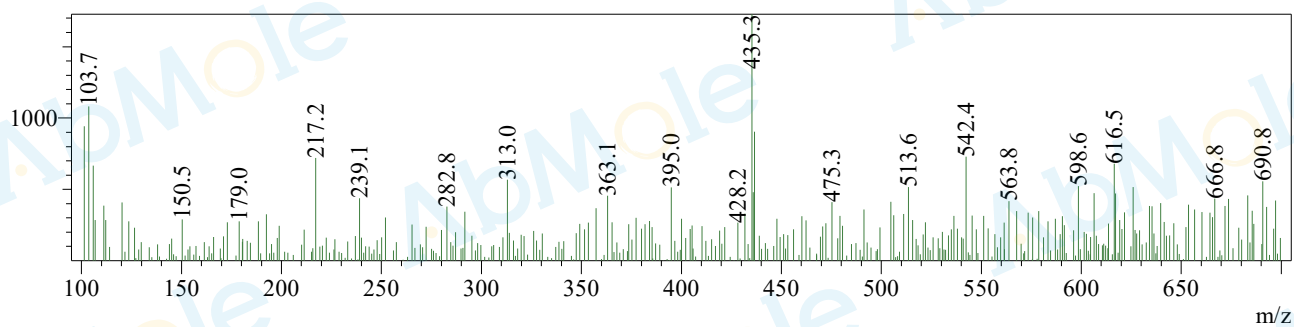


Rt=0.979min 正

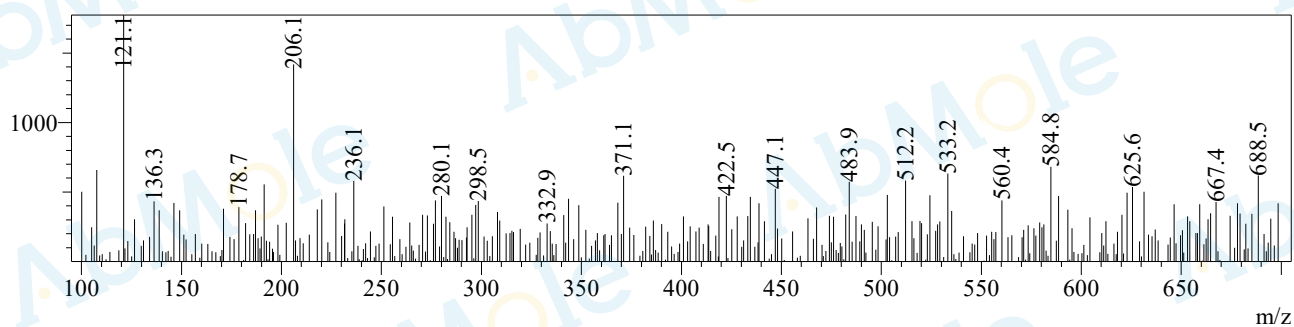

Rt=1.182min 正

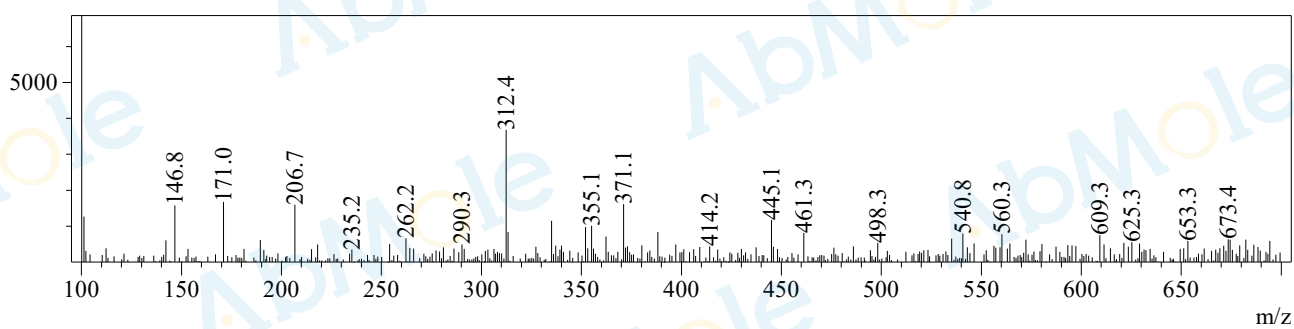

Rt=1.322min 负

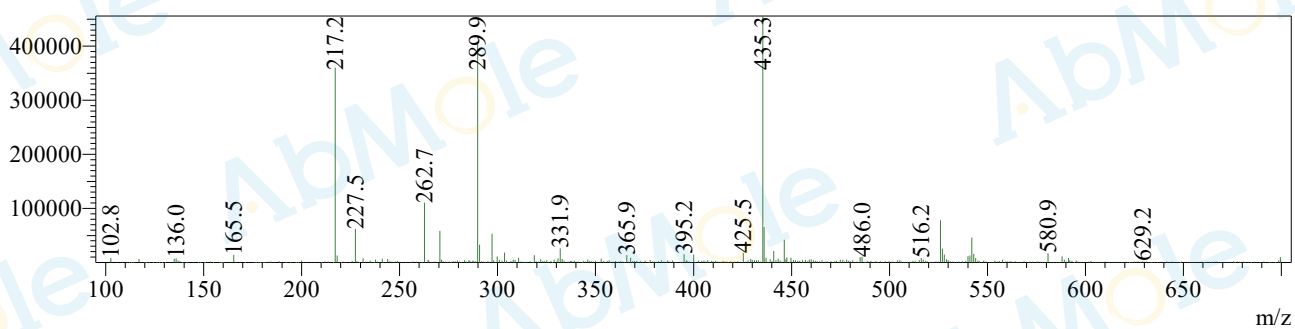

Rt=1.317min 正

