



Sample Name : M5838-Oseltamivir acid-02-D20
Acq. Operator : C22
Spec. Reported : MS Integration
Acq. Method : D:\METHODS\1-POS-3MIN.M
Analysis Method : D:\METHODS\1-POS-3MIN.M
Sample Info : Easy-Access Method: '1-POS-3MIN.M'
Method Info : Column:Xbridge C18(2) (4.6x 50 mm, 3.5um)
Mobile phase: H2O(10 mmol NH4HCO3) (A) / ACN(B)
Elution program: Gradient from 5 to 95% of B in1.5min at 1.8ml/min
Temperature: 50 °C
Detection: UV (214 , 4 nm) and MS (ESI, Pos mode ,103 to 800 amu)


*DAD1 A, Sig=214,4 Ref=off

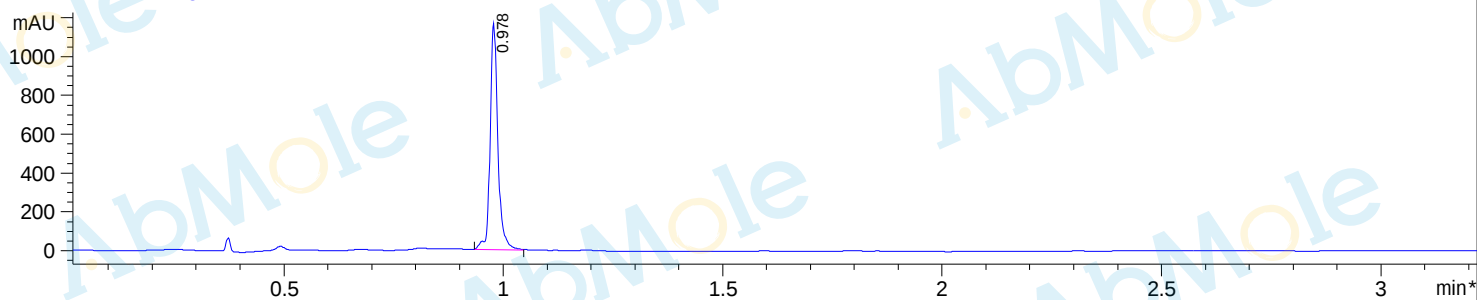

*DAD1 B, Sig=254,4 Ref=off

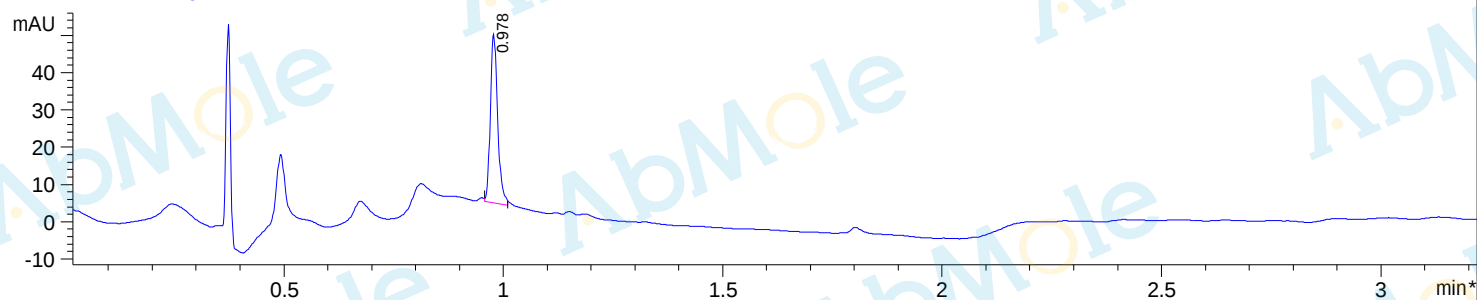

MSD1 TIC, MS File ES-API, Pos, Scan, Frag: 80

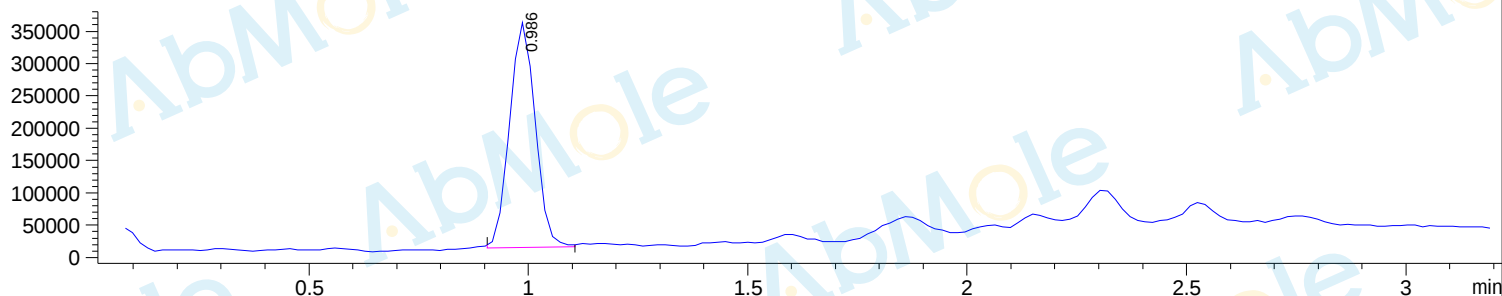

Integration Results for DAD1 A, Sig=214,4 Ref=off

RetTim	Width	Area	Height	Area%
0.98	0.02	1447.77	1164.75	100.00

Integration Results for DAD1 B, Sig=254,4 Ref=off

RetTim	Width	Area	Height	Area%
0.98	0.02	52.76	45.28	100.00

Ret. Time: 0.99

<<< POSITIVE SPECTRA >>>

*ES-API Positive

