

M7790-Emodin-01-LCMS

Sample Information

LCMS


Acquired by : System Administrator
Sample Name : M7790-Emodin-01
Method File : B 5%-95%-P-new - 2.6min.lcm

Date Acquired : 2021/7/27 11:00:01
Injection Volume : 2
Vial : 40


Column : Waters SunFire C18 50*4.6mm 5um 2.000ml/min 2.6min
Pump A : 0.1% FA and 10% ACN in H2O
Pump B : 0.1% FA and 10% H2O in ACN

Chromatogram


mAU

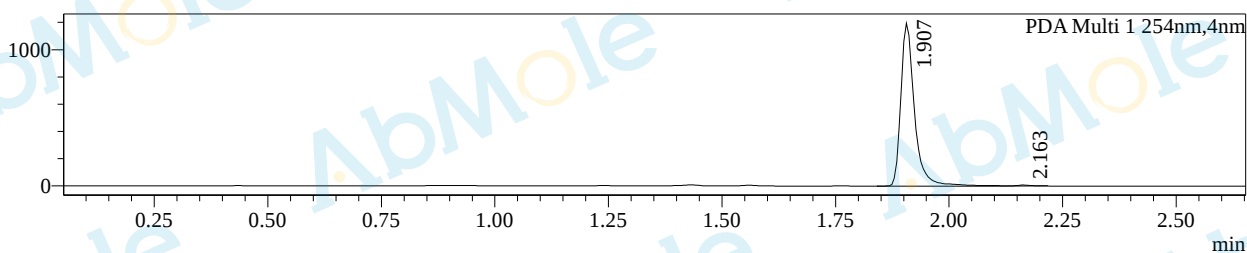

mAU

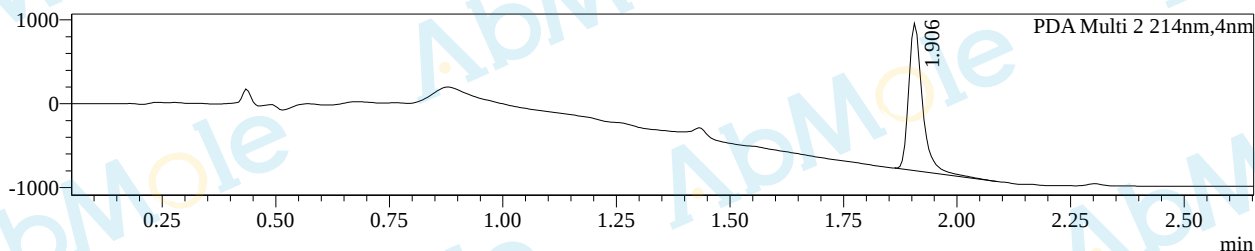

段号1

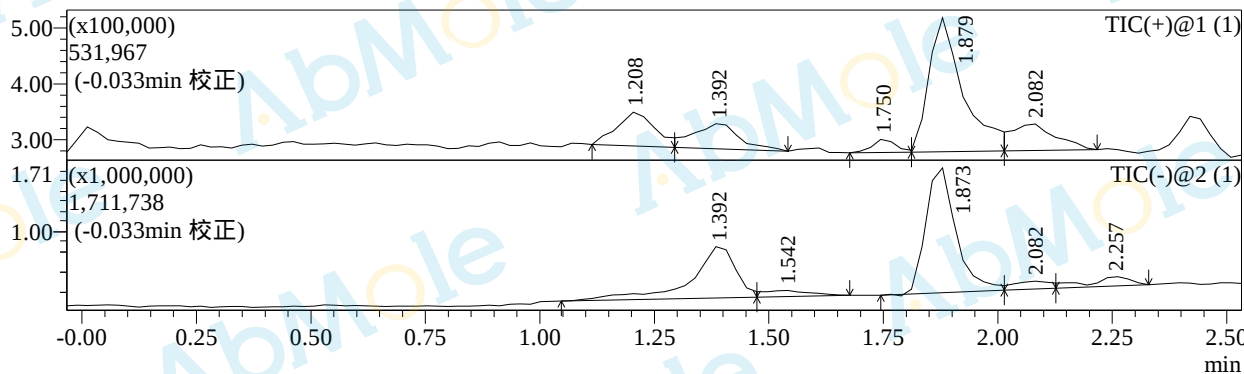

Integration Results for DAD

PDA Ch1 254nm


Peak#	Ret. Time	Height	Area	Area%
1	1.907	1196818	2490513	99.262
2	2.163	7372	18510	0.738
总计		1204189	2509022	100.000

PDA Ch2 214nm


Peak#	Ret. Time	Height	Area	Area%
1	1.906	1758123	3703515	100.000
总计		1758123	3703515	100.000

Mass Spectrum


Rt=1.204min 正

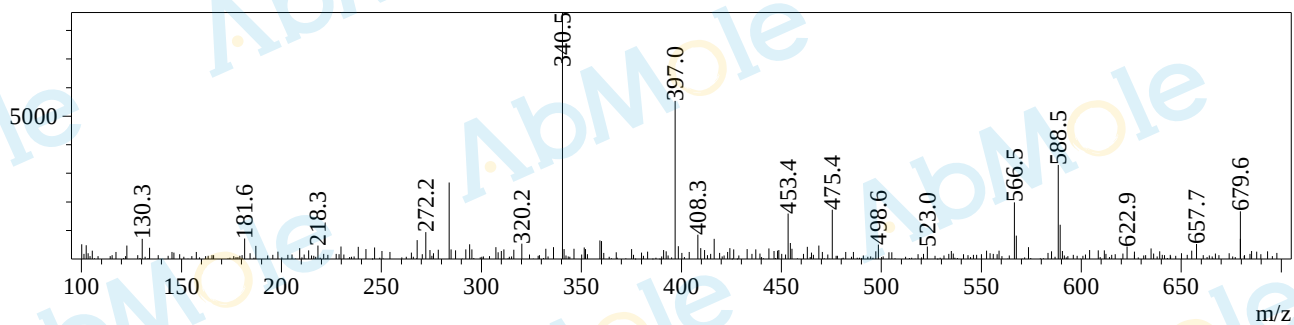


Rt=1.390min 负

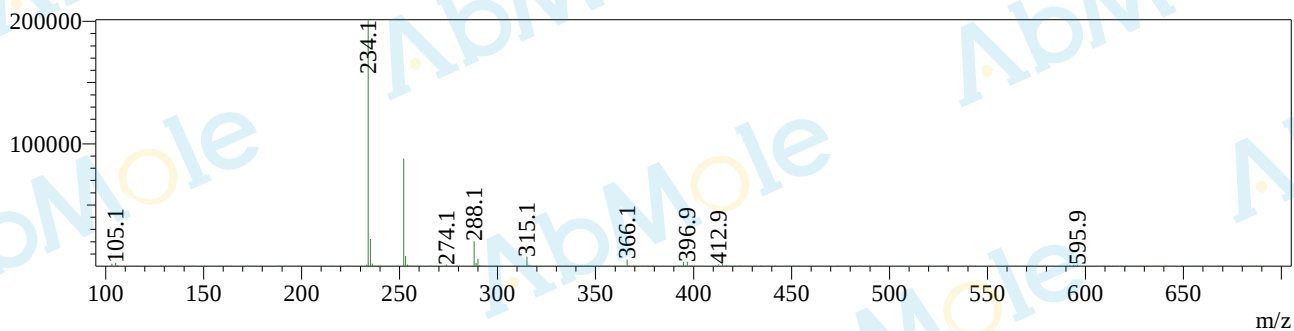


Rt=1.384min 正

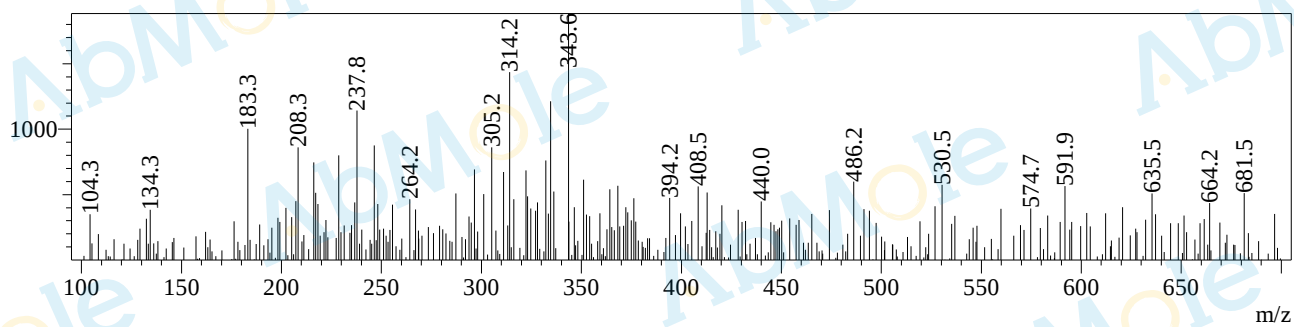


Rt=1.547min 负

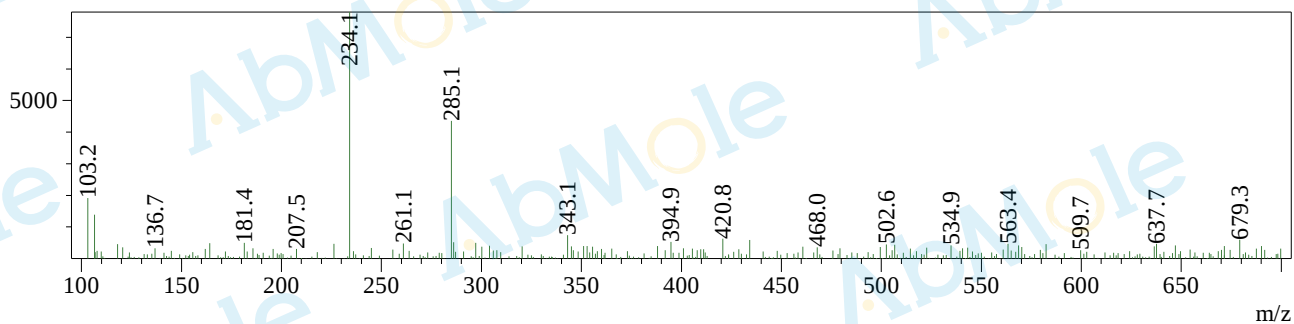


Rt=1.744min 正

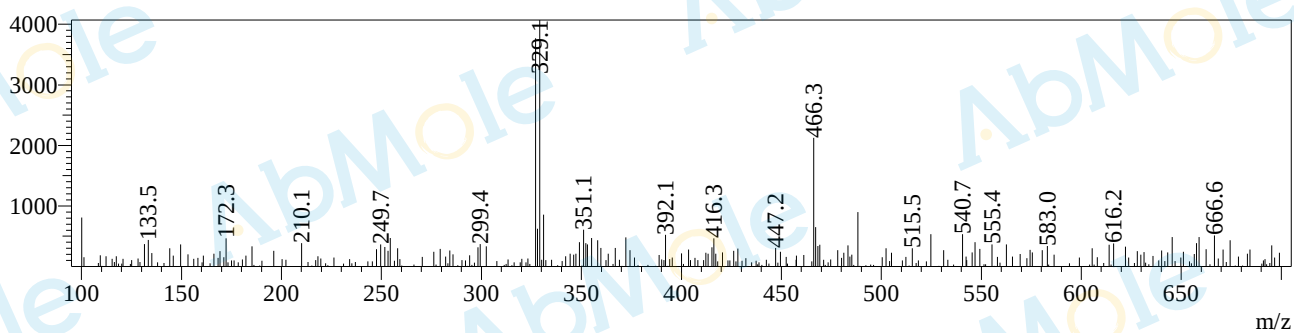


Rt=1.862min 负

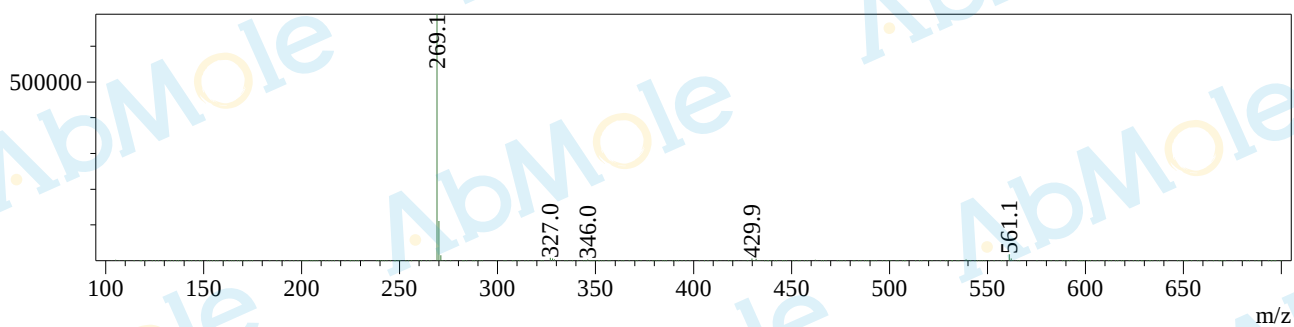


Rt=1.857min 正

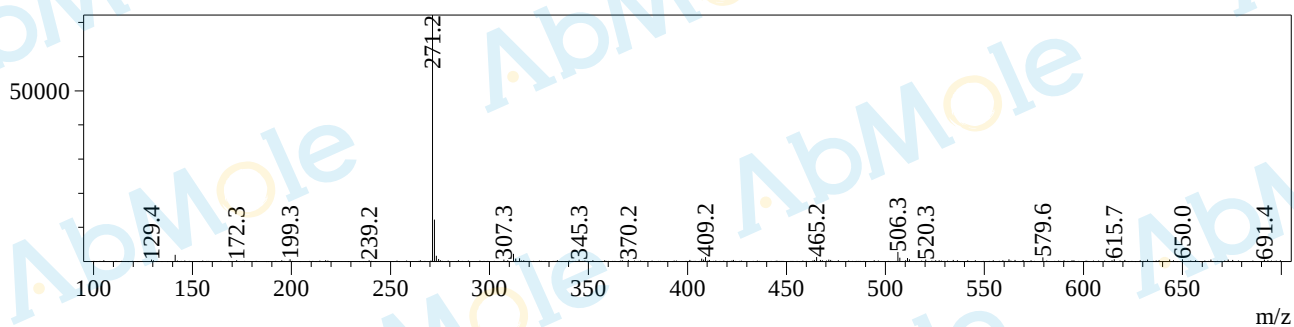

Rt=2.082min 正

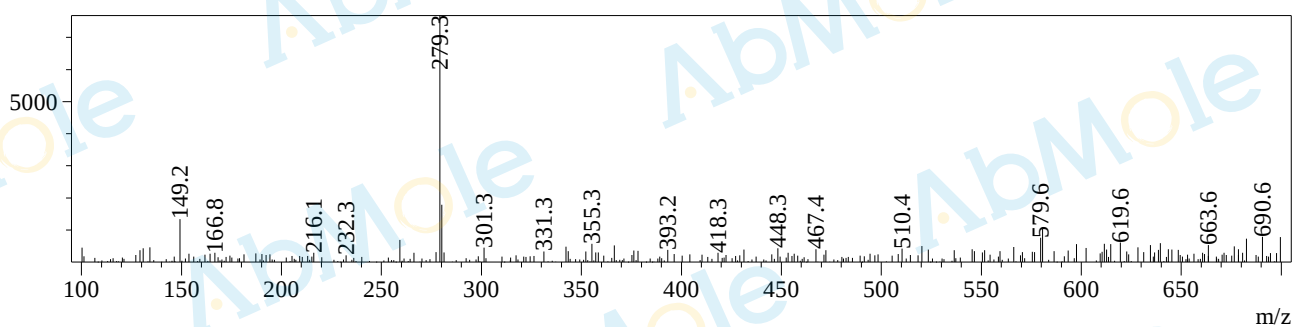

Rt=2.087min 负

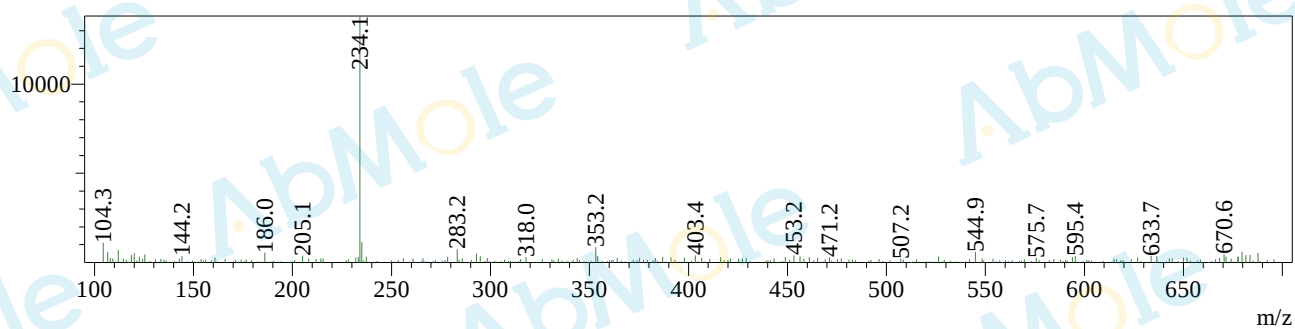

Rt=2.245min 负

