

Shimadzu LCMS-029 Report

Sample Information

LCMS


Acquired by : System Administrator
Sample Name : M9689-Seclidemstat-01
Method File : B 5%-95%-P-new - 2.6min.lcm

Date Acquired : 2021/8/21 10:26:42
Injection Volume : 1
Vial : 71


Column : Waters SunFire C18 50*4.6mm 5um 2.000ml/min 2.6min
Pump A : 0.1% FA and 10% ACN in H2O
Pump B : 0.1% FA and 10% H2O in ACN

Chromatogram


mAU

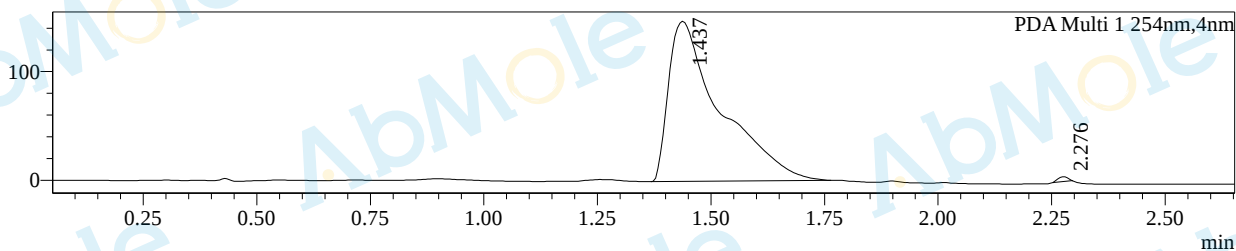

mAU

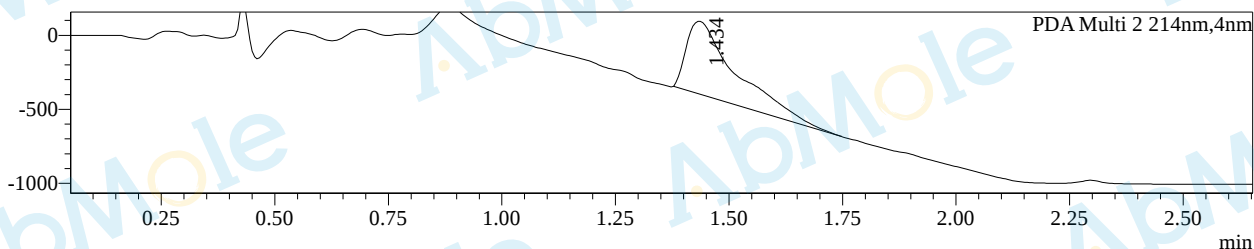

段号1

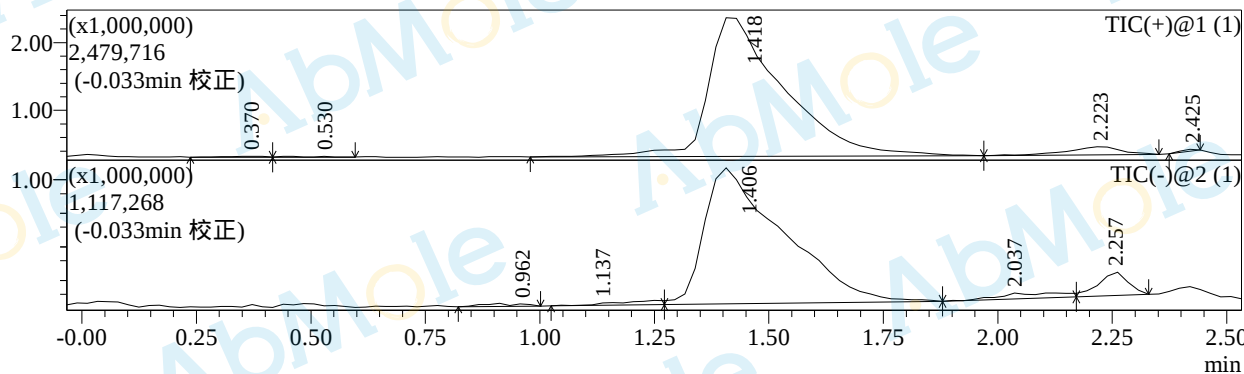

Integration Results for DAD

PDA Ch1 254nm


Peak#	Ret. Time	Height	Area	Area%
1	1.437	147368	1169112	99.420
2	2.276	4371	6825	0.580
总计		151740	1175937	100.000

PDA Ch2 214nm


Peak#	Ret. Time	Height	Area	Area%
1	1.434	490885	3722242	100.000
总计		490885	3722242	100.000

Mass Spectrum


Rt=0.349min 正

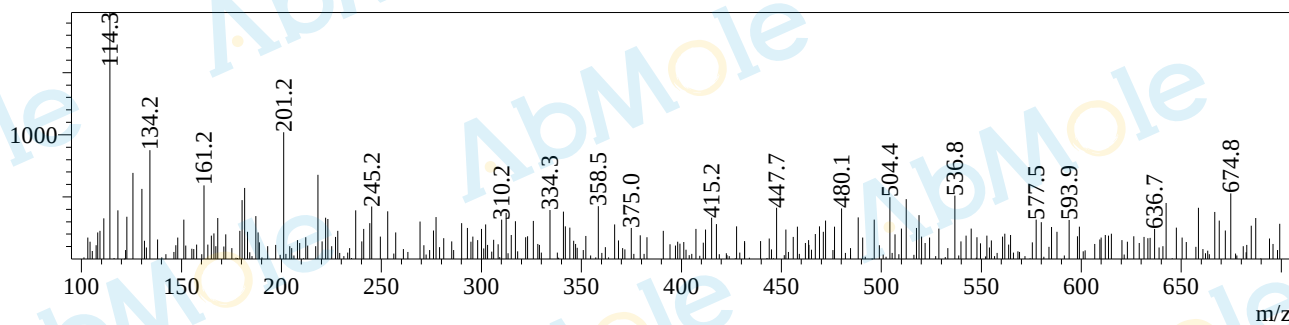


Rt=0.529min 正

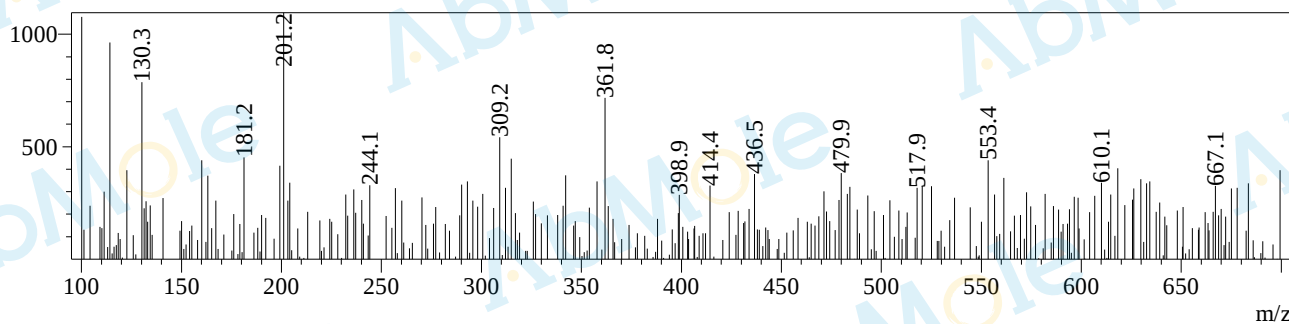


Rt=0.962min 负

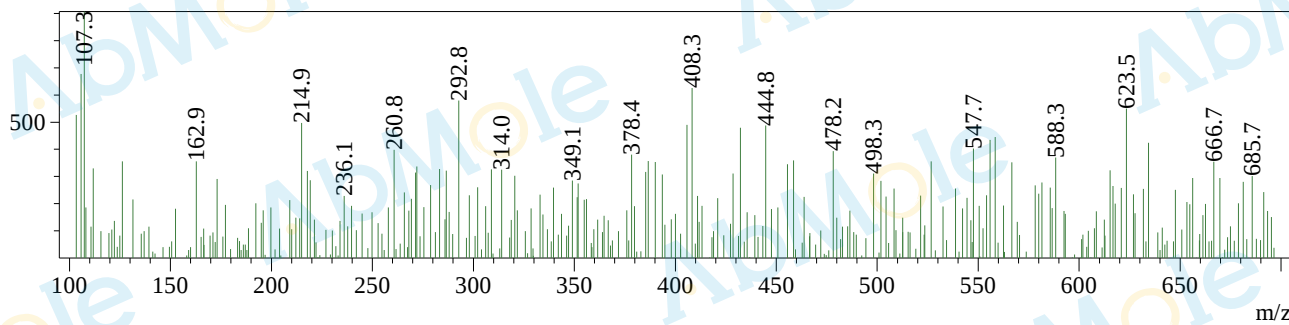


Rt=1.142min 负

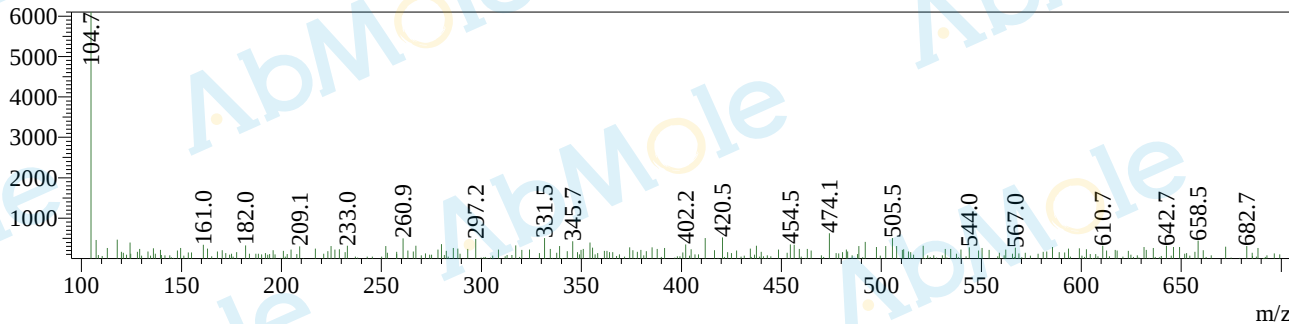


Rt=1.390min 负

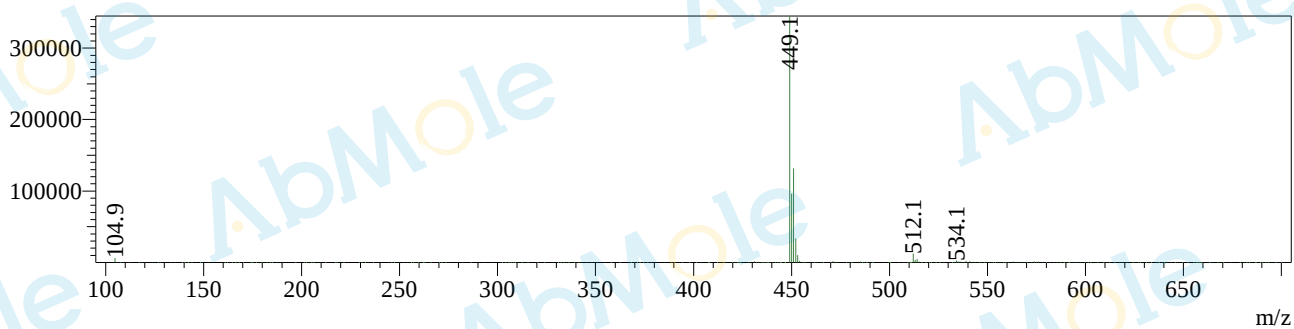


Rt=1.407min 正

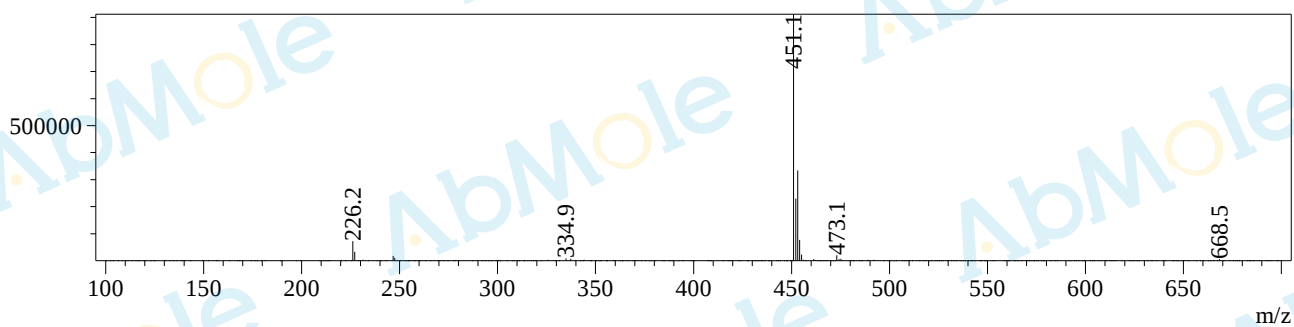


Rt=2.042min 负

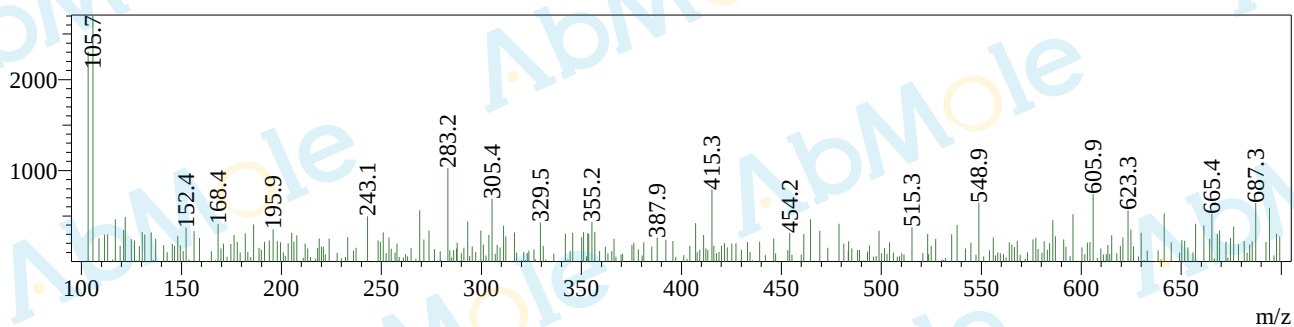

Rt=2.217min 正

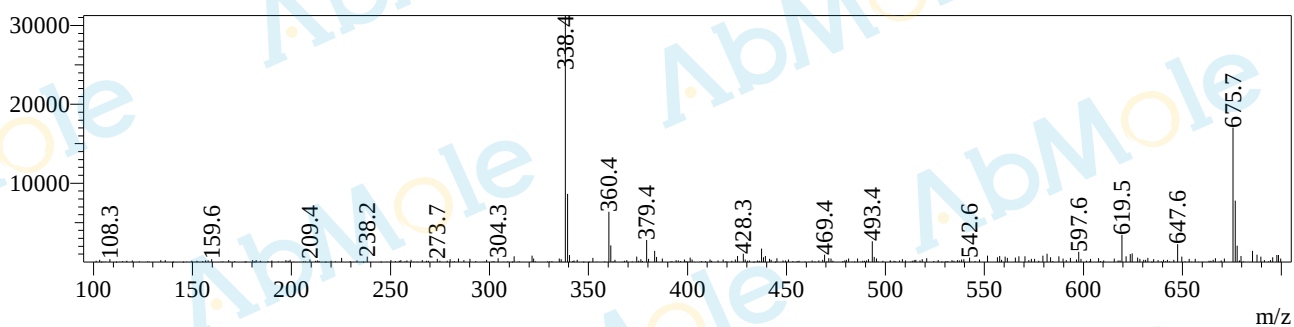

Rt=2.245min 负

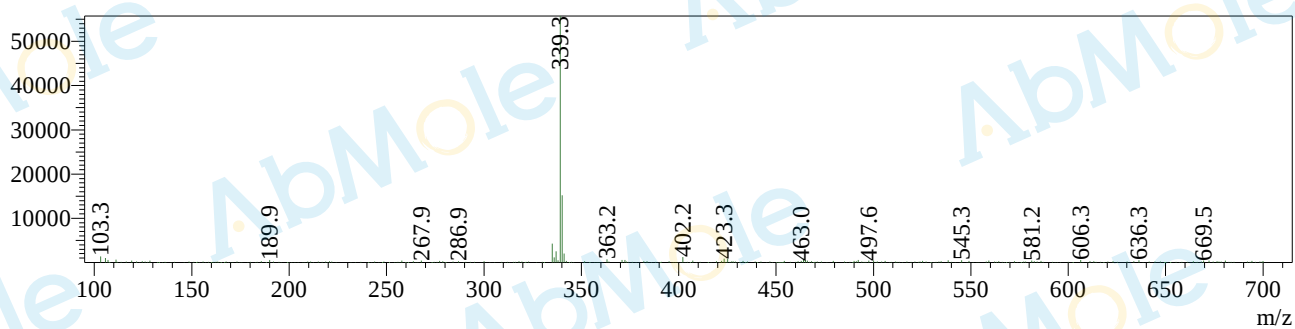

Rt=2.419min 正

