

File D:\DATA\Mitoxantrone dihydrochloride-HPLC

Acq. Method : D: HPLC

Analysis Method : D: HPLC

Sample Info : Easy-Access Method:

Method Info : Mobile Phase: A: water(0.01%TFA) B: ACN (0.01%TFA)

Gradient

Time 0 8 15


B % 5 95 95

Flow Rate: 1.0mL/min


Oven temperature: 40C

Column: SUNFIRE C18, 150*4.6mm, 3.5um

DAD1 A, Sig=214,4 Ref=off

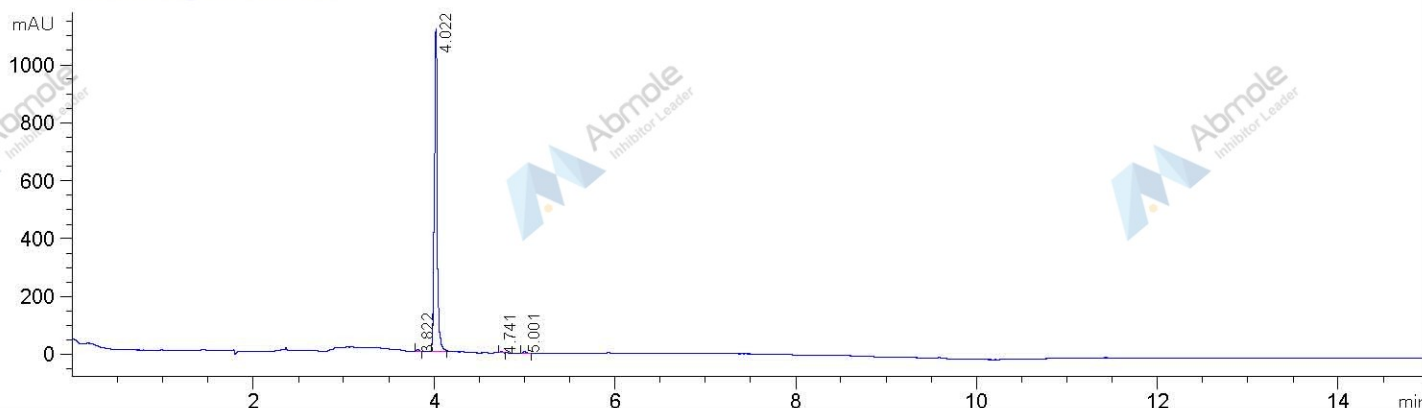

DAD1 B, Sig=254,4 Ref=off

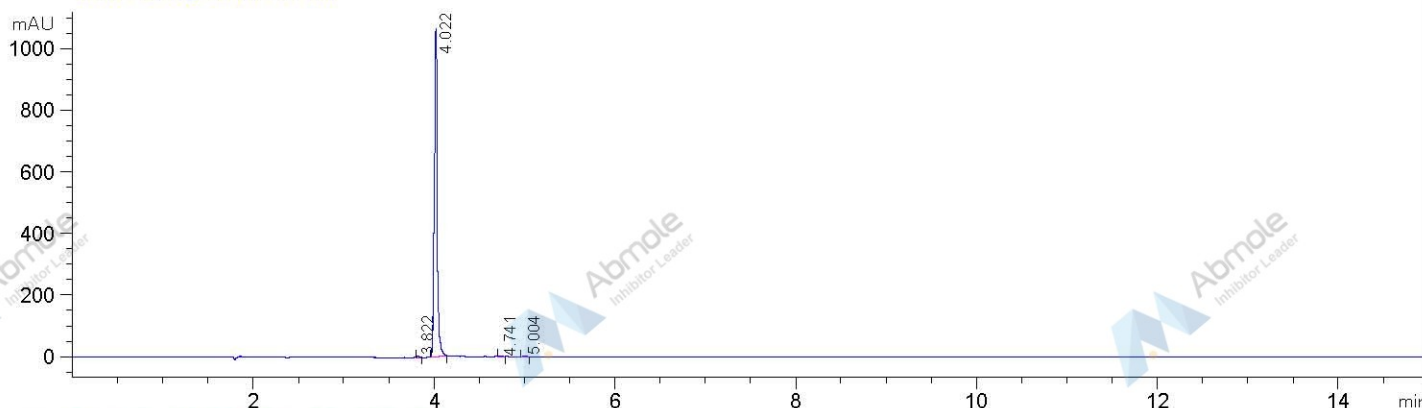

Integration Results for DAD1 A, Sig=214,4 Ref=off

RetTim	Width	Area	Height	Area%
3.82	0.03	11.31	6.50	0.45
4.02	0.03	2472.28	1113.57	98.58
4.74	0.03	7.93	4.29	0.32
5.00	0.04	16.28	6.32	0.65

Integration Results for DAD1 B, Sig=254,4 Ref=off

RetTim	Width	Area	Height	Area%
3.82	0.02	9.10	5.48	0.38
4.02	0.03	2376.92	1062.85	98.84
4.74	0.03	6.55	3.95	0.27
5.00	0.04	12.21	4.83	0.51