

File D:\DATA

Injection Date : 21 Oct 14 3:45

Sample Name : Bevirimat-01

Acq. Operator : XXXX

Acq. Method : C:\Chem32\1\METHODS\1-POS-2MIN.M

Analysis Method : C:\Chem32\1\METHODS\1-POS-2MIN.M

Sample Info : Easy-Access Method: '1-POS-2MIN'

Method Info : Mobile Phase: A: Water(0.01%TFA) B: ACN(0.01%TFA)

Gradient: 5%B increase to 100%B within 0.8min,100%B for 1.1min.


Flow Rate:1.5mL/min

Column: Chromolith Fast Gradien RP-18e,50-3mm.


Column Temperature: 40C

Mass Range:100-1100


*ADC1 A, ELSD

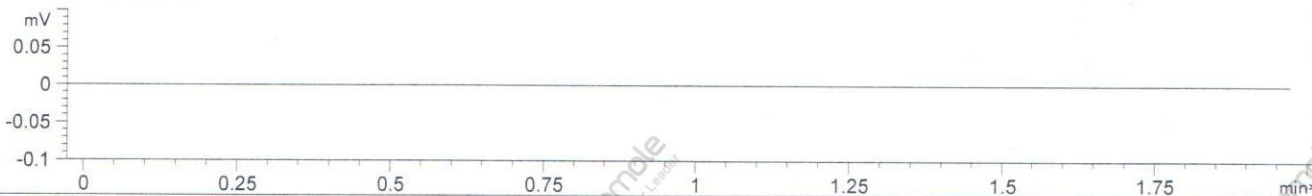


*DAD1 A, Sig=214,4 Ref=off

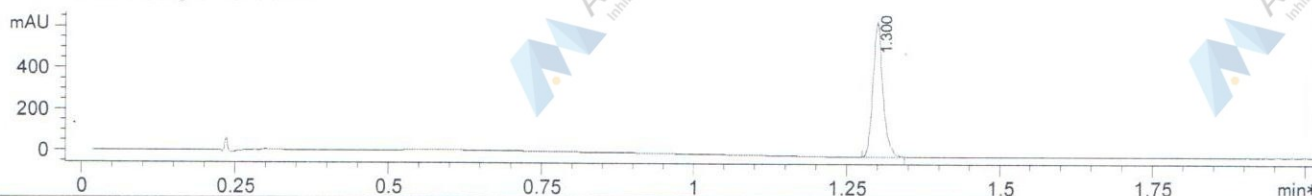

*DAD1 B, Sig=254,4 Ref=off

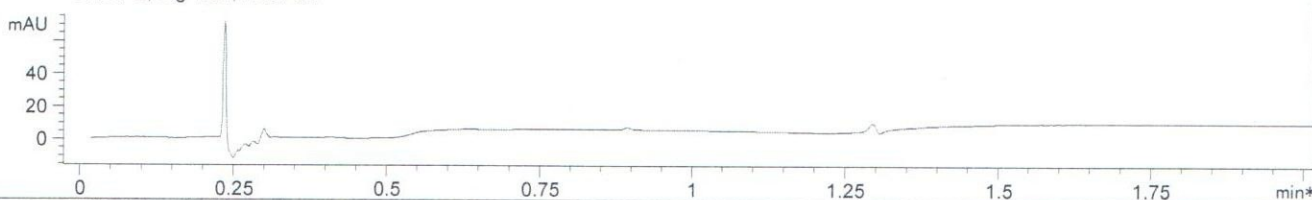

MSD1 TIC, MS File ES-API, Pos, Scan, Frag: 70,

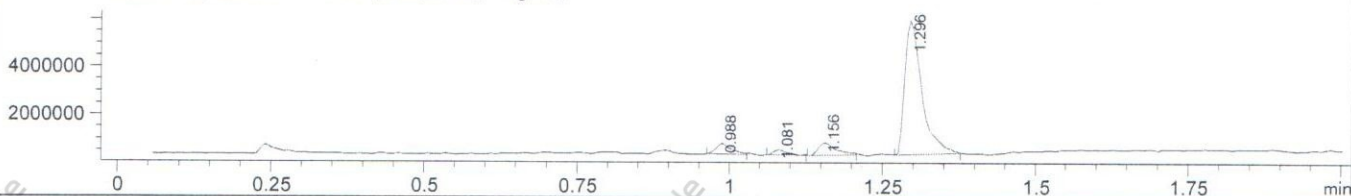

Integration Results for ADC1 A, ELSD

RetTim	Width	Area	Height	Area%
1.30	0.02	758.00	657.55	100.00

Integration Results for DAD1 A, Sig=214,4 Ref=off


RetTim	Width	Area	Height	Area%
1.30	0.02	758.00	657.55	100.00

Integration Results for DAD1 B, Sig=254,4 Ref=off


RetTim	Width	Area	Height	Area%
0.25	0.02	758.00	657.55	100.00

Ret. Time: 0.99


<<<< POSITIVE SPECTRA >>>>

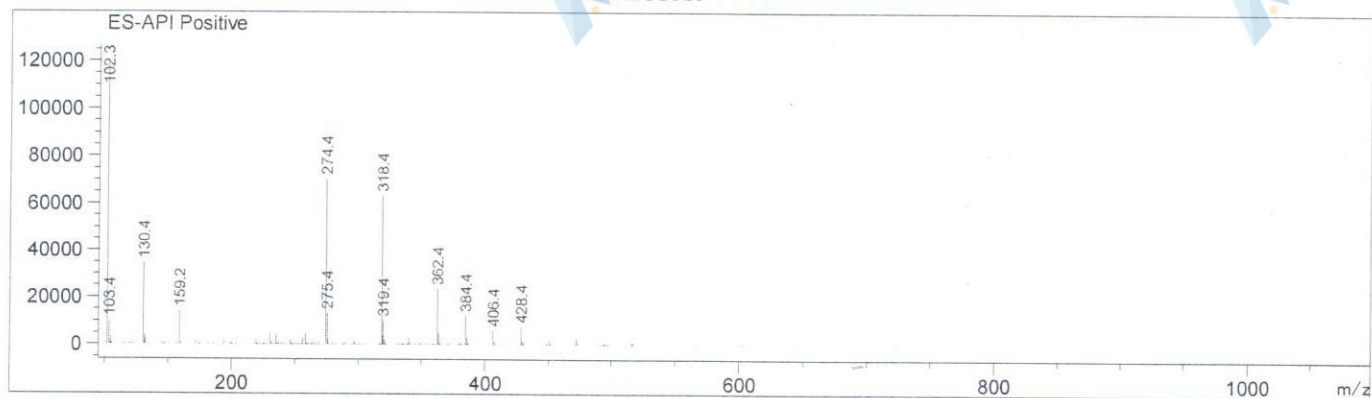

Ret. Time: 1.08

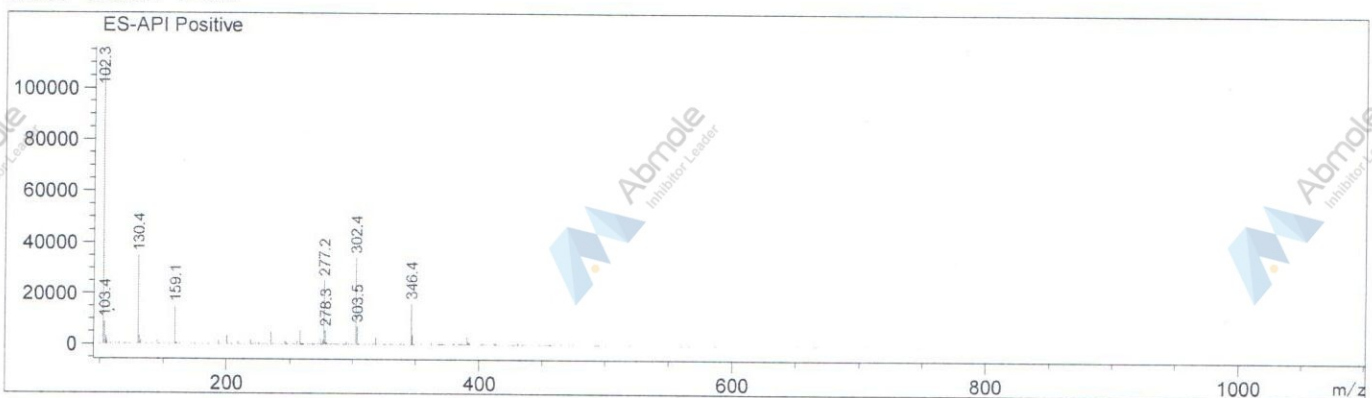

Ret. Time: 1.16

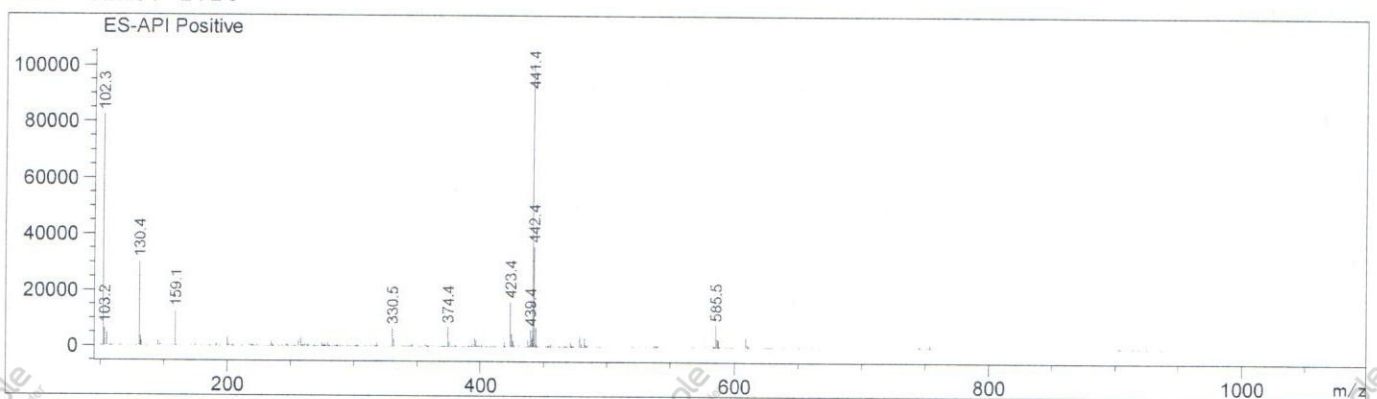

Ret. Time: 1.30

